

Saudi Mais Co. for Medical Products Work Instruction

Title
Production & Test Procedure for Semi Auto Biopsy Needle

Doc. #	:	WI-8-2-2-53
Rev. #	:	0
Date. #	:	17.10.2020
Page. #	:	1 / 4

1. Purpose:

To stipulate the required tests for Semi Auto Biopsy needle to ensure good quality products before releasing to the market.

2. Scope:

This procedure is applicable for product manufacturing process of Semi Auto Biopsy needle by Saudi mais Co, for medical products.

3. Definition:

Semi Auto Biopsy Needle: A medical device which is *Semi-automatic* spring loaded *biopsy* system designed to cover a wide range of soft tissue core *biopsy* procedures. Lightweight compact instrument specially designed to perform *biopsies* under CT scan guidance.

4. Responsibilities:

- 4.1 Quality Controller is responsible for sampling & testing.
- 4.2 Quality control supervisor is responsible for checking the obtained results.
- 4.3 Machine operators and assembly peoples are responsible to produced good quantity and good quality products.

5 Procedure:

5.1 Receiving Bulk components :

- Ensure the receiving slip from supplier is correct gauge of biopsy needle as per our P.O.
- Inspect the carton condition, must be in good condition, no damaged at all cartons with products.
- Inspect the products condition & components which written on receiving slip from supplier.
- Take samples & inspect in order to AQL procedure, product description, specification with apply the AQL sampling plans for normal inspection.
- **Needle Protector**
Quality controller is responsible to check the needle protector it must be perfectly fit with needle and cover the needle to protect it.
- **Needle**
Quality controller is responsible to check the needle it must be stainless steel, smooth surface, sharp tip needle, no tip damage, no blockage and free from irregularities defects. Record the results into OBM inspection report.
- **Spring**

Prepared by : Q.S
Esmayaden Alango

Approved by Q.C.M:
Mr. Javeed Gafoor

Controlled Copy

Audi Mais Co. for Medical Products

Work Instruction

Title
Production & Test Procedure for Semi Auto Biopsy Needle

Doc. #	:	WI-8-2-2-53
Rev. #	:	0
Date.#	:	17.10.2020
Page.#	:	2 / 4

Quality controller is responsible to check the spring it must be flexible not too hard not too loose.

- **Graduation Scale**

Quality controller is responsible to check the graduation marks on the needle for accurate measurement, 1 cm distance in every mark.

- **Safety lock**

Quality controller is responsible to check the safety lock must be good function.

- **Hand Grip**

Quality controller is responsible to check the hand grip of the product that it can be easily to grip & good function.

6. **Components Assembly**

- Production controller is responsible to issue required gauge biopsy needle, stylet needle, needle holder, stylet holder, and protectors.
- Production controller is responsible to issue production order for Assembly of biopsy needle.
- Quality controller is responsible to check the gluing in stylet holder and check by pull test and make sure that stylet is strong not pullout.
- Quality controller is responsible to check the assembled products by randomly and make sure all are proper inserted into connector, no excess glue, & should be strong gluing, no pull out and record the results into In- process inspection report.
- Quality controller is responsible to perform the following tests in order to tensile test for needles and connectors gluing.

6.1 **Checking for needle gluing (Tensile test) by Zwick Machine.**

- Insert the needle into apparatus gripper of universal testing machine.
- The test will be performed as per table requirements.

Needle OD (mm)	Min. force of gluing
0.6 - 0.9 mm	25 N
0.9 - 1.2 mm	30 N

Prepared by : Q.S
Esmayaden Alango

Approved by Q.C.M:
Mr. Javeed Gafoor

Controlled Copy

audi Mais Co. for Medical Products

Work Instruction

title

Production & Test Procedure for Semi Auto Biopsy Needle

Doc. #	:	WI-8-2-2-53
Rev. #	:	0
Date.#	:	17.10.2020
Page.#	:	3 / 4

1.2 - 2.8 mm

40 N

6.2 Visual Inspection.

- The needle must be 100% stainless steel material and smooth surface.
- Check the needle orientation and bevel.
- Check the needle holder are tightly assemble into needle which is comply with standard BS EN ISO 10555-5 : 2013
- Check the needle tip end and must be sharp and free from irregularities defect, no edge, no tip bend.
- Check the needle dimension as per our required specification and size.
- Check the needle eye, Quincke bevel type.
- Checked the Needle should be free from any blockage.

6. Packaging and Labeling:

- Production controller is responsible to issue production order for packaging.
- Machine operator is responsible to fix die with clichés and set approved parameters for packaging.
- Each product package in blister pack, the label contains the following.
Product name and code no., reference no., size, lot no., type, temperature, logo, sterilization logos, production date, expiry date, etc.
- Machine operator is responsible to give initial sample for inspection and quality controller is responsible to check the package it must be strong sealing.
- Then to pack in specified small boxes, and arrange in big carton.
- The big carton must be with label and detail information.
- During packaging process, quality controller to be check the package sealing, peeling, sealing width, forming, cutting and printing.


7. Sterilization:

- The finished products are sterilized as per work instruction (WI-7-5-1-1).


8. Final Inspection and Product Release

- The Quality controller must be check, the unloaded products after sterilization.
- The quality assurance manager or quality control manager have to review the Completed cycle parameters and inspection reports both documents have to be sign.

Prepared by : Q.S
Esmayaden Alango

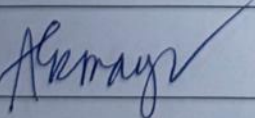

Approved by Q.C.M:
Mr. Javeed Gafoor

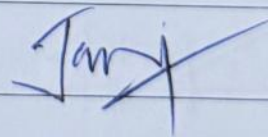

Controlled Copy

audi Mais Co. for Medical Products

Work Instruction

title

Production & Test Procedure for Semi Auto Biopsy Needle

Doc. # : WI-8-2-2-53

Rev. # : 0

Date.# : 17.10.2020

Page.# : 4 / 4

- After confirmation all the test from the Lab, to release the sterilized products, the document once again have to review the quality control manager then sign.
- The production department responsible to record the quantity of products before transferring to the finished goods ware-house and then it could be released to the customer.

9. Sampling:

- Inspection and test report for Raw Material : F-823 - T90
- Inspection & Test report for needle and Stylet : F-823 - T58
- In-process control for packaging machine and Final repacking : F-823-i4
- Final inspection and Release Record : F-823 - F23


10. Reference:

- BS EN ISO 9626:2016 Stainless steel needle tubing for the manufacture of medical devices.


11. Modification / Documentation change History:

Date	Changes
17.10.2020	1 st Issue

Prepared by : Q.S
Esmayaden Alango

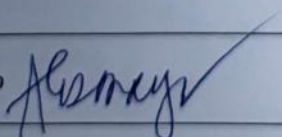

Approved by Q.C.M:
Mr. Javeed Gafoor

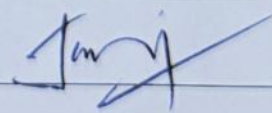

Controlled Copy